RGMS-30

QUICK OPERATION GUIDE

*Please refer to the IFU for the contraindications or precautions.

Components



(1) Transmitter (2) Auxiliary Tool (3) Disposable Glucose Sensor (4) App

❶ Prepare for the Applying

(1) Download the App

For Android:

Search the app "MQ CGM" in Google Play or you can scan the QR code to download the app.

For iOS: Search the app "MQ CGM" in App Store.

Register and login according to the app steps.

(2) Click "Start monitoring"



Note : please turn on the permissions of Bluetooth, Location, Camera, Photos and Videos, Music and Audio, and Nearby Devices during the session. Please refer to the "Help" in "my" menu to ensure the realtime performance of the app.

(3) Scan the Code of the Sensor



Scan the code

Select "Scan code" or "Input" in app to identify the ID of the sensor.

❷ Select and Clean the Sensor Site

Clean the Sensor Site with Alcohol Wipe




Apply to upper arm (avoid muscle). Clean skin with alcohol swab before use.

❸ Assemble the CGM



❹ Connect with App

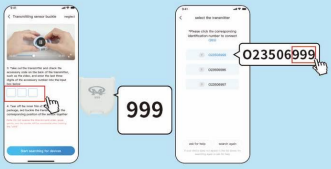

1.Enter the 3-digit code on the transmitter.

2.Search in app and select the correct device to connect.

Note: Please do not wear the device if it could not be connected.

❺ Asseble CGM with Auxiliary Tool



push down to the end | pull out | unlock the auxiliary tool clockwise

❻ Wearing the Device



Tighten the skin of the sensor site and place the auxiliary tool closely over the site. | Press the white button | Keep it still for 5 seconds and move it

Press "Applying completed" in app. You will see your glucose reading there after warming up of 1 hour. Note: The CGM is factory calibrated. If you found significant difference from your BGM reading, calibrate it when your glucose level is stable. (best in morning on empty stomach).

Remote Monitor



Tips: VIP LINE 400–0572–278

1. Please keep the mobile bluetooth function on during the wearing.
2. Please use the same fingertip blood glucose meter in the sensor session.
3. Please do not press End Monitoring in app during monitoring.
4. The transmitter could be used repeatedly for 3 years. Please keep it well after separation with sensor.

CE-3OZN-02V.01

End Monitoring

